REDLIGHT

SVĚTELNÁ TERAPIE – MINI PANEL

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BEZDRÁTOVÉ

KOMPAKTNÍ

**BATERIE
5 200 mAh**

144 LED DIOD

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144 LEDS

CONTENTS

1. Warnings and Disclaimer	55
2. How Red Light Therapy Works	57
3. Wavelength Overview	58
4. Package Contents	59
5. Control Overview	60
6. How to Use	61
7. Product Specifications	63
8. Notes	63
9. Disclaimer	64
10. EU Compliance	64
11. Safety Instructions	65

Do not place this product in environments containing corrosive chemicals, as corrosion may cause product failure.

1. WARNINGS AND DISCLAIMER


- Before use, we recommend consulting your healthcare provider if you believe you are sensitive to light. Certain cold and allergy medications, pain relievers, and medications used to treat infections may cause light sensitivity.
- Before use, we recommend consulting your healthcare provider if you are pregnant, have any suspicious or cancerous lesions, have recently received a steroid injection, or are using topical steroids.
- If you experience discomfort or have any concerns regarding the light, discontinue use immediately and contact your healthcare provider.
- Disconnect the device from the power source when it is not in use.
- This device is not intended to diagnose, treat, cure, or prevent any disease. The company makes no express or implied warranties regarding the therapeutic effectiveness of this product in relation to physical health, skin conditions, mental health, or any other medical conditions. Any use of this product for therapeutic purposes must be carried out under the strict supervision of qualified healthcare professionals and in accordance with their prescribed instructions. Users are advised to fully understand the intended scope of use of the product. Any legal liability arising from improper operation, related health risks, or subsequent claims shall be borne solely by the responsible party.

To maintain the service life and functionality of the device, always store it in a cool, dry place. For proper use and safety, please read this manual carefully and follow all instructions.

Do not do any of the following:

- Do not allow children to use the device without adult supervision.
- Do not use the device on infants.
- Do not allow sleeping, disabled, or unconscious persons to use the device.
- Do not exceed the maximum recommended usage time of 30 minutes.
- Do not use on open or fresh wounds.
- Do not use in combination with ointments or balms containing heat-producing ingredients. Doing so may result in skin burns.
- Do not modify or attempt to repair this device.
There are no user-serviceable parts inside.
- Do not clean electrical components with water or other liquids.
- Ensure that the device and its components, such as the battery, do not fall into water or any other liquid.
- Do not place this product in environments with excessively high temperatures. High temperatures may cause damage to or explosion of the battery.
- Do not store or operate the device at temperatures above 48°C (118°F).
- Do not use the device if it is damaged. Continued use of a damaged device may result in electric shock or injury.
- Do not charge the battery while using the product.
- During use, protect the device from drops, impacts, or other mechanical damage.
- If repairs are required, please contact an authorized service center or your dealer.
- This product contains precision electronic circuits. Unauthorized disassembly may result in malfunction.
- It is recommended to use original cables only. Any damage caused by the use of non-original accessories is not covered by the warranty.


2. HOW RED LIGHT THERAPY WORKS

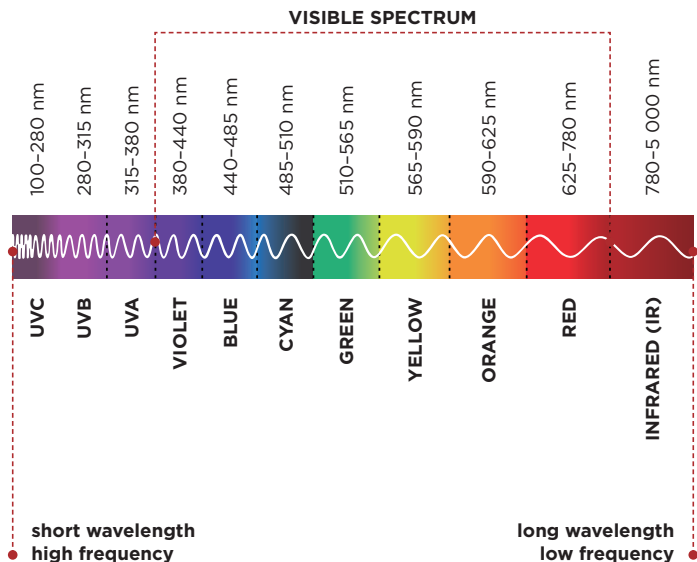Specific photons found in red and near-infrared light interact with a photoreceptor in our cells called cytochrome c oxidase. This interaction helps the mitochondria in our cells use oxygen more efficiently, allowing them to produce more ATP (adenosine triphosphate). Adenosine triphosphate is a complex organic chemical that provides energy to power many processes in living cells, such as muscle contraction, nerve impulse transmission, and chemical synthesis. ATP, which is found in all forms of life, is often referred to as the “molecular unit of currency” of intracellular energy transfer. When our cells have more energy, they simply function better, and as a result, the body functions better as well.


3. WAVELENGTH OVERVIEW

- **Red Light:** Visible to the human eye, red light has wavelengths between 600–700 nanometers (nm) and is most effective at the surface of the skin. Red light may help improve skin health, support collagen production, and promote hair growth.
- **Near-Infrared Light:** Invisible to the human eye, near-infrared light has wavelengths between 700–1200 nm and can penetrate deeper into the body, affecting muscles, joints, and bones.

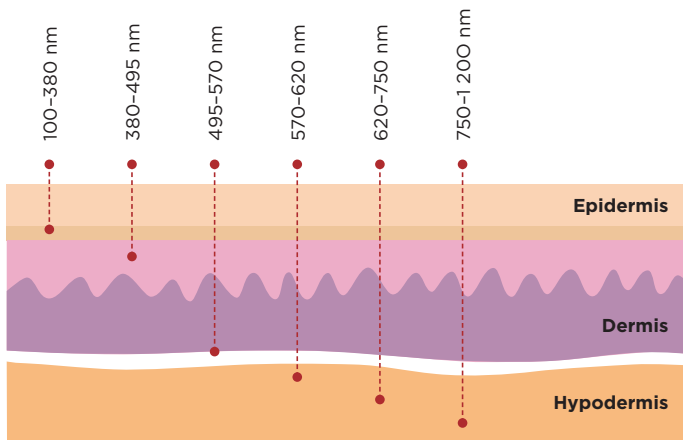Wavelength	Color Spectrum
100–380	UV
380–495	Blue Light
495–570	Green Light
570–620	Yellow Light
620–750	Red Light
750–1200	NIR

4. PACKAGE CONTENTS

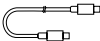Accessory	Quantity	Image
Mini Panel	1	
Cover / Shade	1	
USB-C Power Cable	1	
Power Adapter	1	
Protective Glasses	1	
Stand	1	
User Manual	1	

Power Adapter

Input: 100-240 V, 50/60 Hz

Output: 12 V / 1.67 A

USB-C Power Cable



Input: 12 V / 1.67 A (max. 20 W); 20 V / 3.25 A (max. 65 W)

Protective Glasses

Do not look directly into the light source (LEDs) during use.

We recommend wearing the protective glasses included in the package during use.



5. CONTROL OVERVIEW

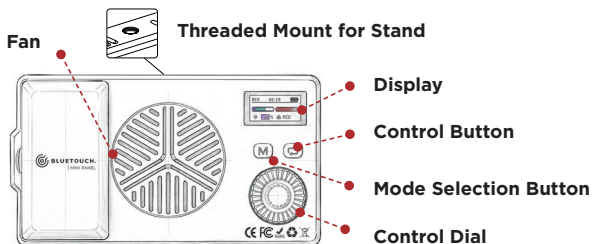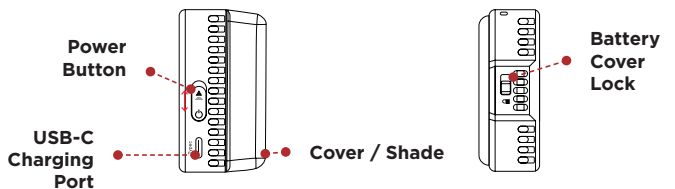6. HOW TO USE


Step 1:

Slide the switch upward to turn on the device. A symbol will appear on the display. Turn the dial clockwise to switch on the light.

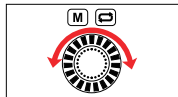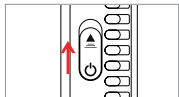Step 2: “NIR/RED Mode”


Briefly press the „“ button to switch between light intensity adjustment and NIR/RED modes. In NIR/RED mode, rotate the dial to increase or decrease the value. (Note: 0% = NIR, 100% = RED, 1-99% = N-R). In Intensity Mode, rotate the dial to increase or decrease the light intensity (1-100%).

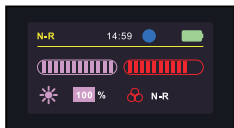Step 3: “Smart Mode”

Briefly press the „“ button to enter Smart Mode. Rotate the control dial to select one of the 5 available modes: (10 Hz, 7 Hz, 0,25 Hz, 5 Hz (Pet Mode), Mute Mode). Briefly press the „“ button to switch to brightness adjustment (1-100%).

Note: The timer cannot be adjusted here. It must be set in the “Settings Mode” (5/10/15 minutes). Mute Mode is fixed at 30 minutes and cannot be adjusted.

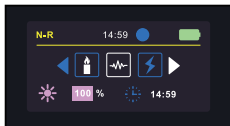10 Hz

5 Hz

7 Hz

20 % (20-30 min)

10 % (10-20 min)

5 % (0-10 min)

0.25 Hz (Pulse Mode)

Note: Long press the „**M**“ button — Entering Mute Mode.

Activates a 30-minute countdown. The fan is dimmed and turned off, and the brightness is set to 20%. The brightness will gradually decrease every 10 minutes until the light turns off completely.


Step 4: “Settings Mode”

Briefly press the „**M**“ button to enter Settings Mode.

The timer can be set to 5, 10, or 15 minutes by rotating the dial.

Briefly press the „“ button to switch to the fan control function.

The fan can be turned on or off.

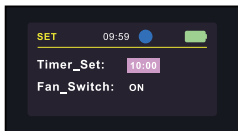Note: Long press the „“ button — Turn the Fan On/Off.

- When the fan is turned off, the maximum brightness is limited to 20%.
- To increase the brightness above 20%, the fan must be turned on.

7. PRODUCT SPECIFICATIONS

Model	BLUETOUCH - Light Therapy; Mini Panel	Charging	USB-C
Wavelength	RED 660 nm / NIR 850 nm	Charging Time	≈ 1.2 h
Charging Voltage	12 V / 1,67 A	Power Supply	Built-in Battery + USB-C Power Supply
Smart Modes	10 Hz, 7 Hz, 0.25 Hz (Pulse Modes), 5 Hz (Pet Mode), Mute Mode	Color	White / Black / Pink
Brightness Adjustment	1-100 %	Weight	174 g
Dimensions	142 × 79 × 30 mm		

8. NOTES

Our company reserves the full right to revise and interpret this manual and makes every effort to ensure the accuracy of its contents. The contents of this document have been carefully reviewed. In the event of any printing errors or omissions, Bluetouch s.r.o. reserves the right of final interpretation. At the same time, Bluetouch s.r.o. reserves the right to improve and modify its products and this manual. Product specifications may change without prior notice. Some images in this manual are for illustrative purposes only. If an image does not correspond to the actual product, please refer to the actual product.

9. DISCLAIMER

1. This product **is not intended** for the diagnosis, treatment, or prevention of any disease. The company makes no claims, statements, or warranties regarding the ability of this product to treat any physical, dermatological, or mental condition. If you have any health issues requiring medical attention, please consult a qualified healthcare professional.
2. The company assumes no legal liability for any claims arising from safety issues caused by **improper** use of the product.
3. **Built-in Rechargeable Battery**
This product contains a built-in rechargeable battery. Do not disassemble, puncture, or expose the battery to fire, water, high temperatures, or mechanical damage. Do not charge the product unattended, on a flammable surface, or near flammable materials. If the product overheats, emits an odor, changes shape, or shows signs of damage, stop using it immediately and contact customer service.

10. EU COMPLIANCE

This product has been designed and manufactured in accordance with the applicable legal regulations of the European Union relating to consumer electronics.

The product complies with the requirements of the following EU directives and regulations:

- Directive 2014/30/EU on electromagnetic compatibility (EMC),
- Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS),
- ČSN EN IEC 55014-1.

Based on compliance with the above requirements, the product is marked with the CE marking. The EU Declaration of Conformity is available from the importer and can be provided upon request to the relevant supervisory authorities or customers. This product is intended exclusively for **wellness and relaxation use**. It is not a medical device within the meaning of Regulation (EU) 2017/745 and is not intended for the diagnosis, treatment, or prevention of diseases.

11. SAFETY INSTRUCTIONS

Safety Instructions

Use the device only for purposes described in the manual. Do not use if the device is visibly damaged, overheated, or emits an unusual smell. Persons with reduced heat sensitivity must use the device only under supervision. Do not use on open wounds, burned skin, fresh scars, infections, or inflamed areas. Do not use if you have an implanted electronic device (e.g., pacemaker) unless approved by a doctor. Keep out of reach of children unless stated otherwise. Do not use while sleeping. Do not expose the device to water or high humidity. Do not use near flammable materials.

Contraindications and Warnings

Not suitable for persons with photosensitive conditions. Not suitable for light-sensitive epilepsy. Not suitable during pregnancy (unless approved by a doctor). If irritation, redness, or pain occurs, stop using the device immediately.

Maintenance and Cleaning

Allow the device to cool down after use. Clean only with a soft dry or slightly damp cloth. Do not use alcohol, solvents, or aggressive cleaning agents. Do not immerse the device in water.

Storage

Store in a dry place away from direct sunlight.
Storage temperature: -10 °C to +45 °C. Protect from dust and humidity.

Device Disposal

Do not dispose of with household waste.
Take to a specialized collection point.

Manufacturer

BlueTouch, s.r.o.

**Zelený pruh 1974/109,
140 00 Prague, Czechia**

ID: 09437223, VAT: CZ09437223

info@bluetouch.cz