REDLIGHT

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REDLIGHT LIGHT THERAPY FOR KNEES

- 3 TIMER MODES
- 6 HEATING LEVELS
- 2 WAVELENGTHS

WIRELESS

**UNIVERSAL
SIZE**

**5 000 mAh
BATTERY**

44 LED LIGHTS

 **BLUETOUCH.**

**WARNING**

This device is not intended for diagnosis, treatment, or prevention of diseases. The company provides no explicit or implied guarantees regarding therapeutic effectiveness related to physical health, dermatological conditions, mental health, or any other medical field. Use of this product for therapeutic purposes must be supervised by qualified healthcare professionals and follow their instructions. Users should fully understand the intended purpose of this product. Any legal responsibility arising from misuse, health risks, or subsequent claims lies solely with the responsible party.

USER MANUAL

BLUETOUCH - LIGHT THERAPY FOR KNEES

Please read this manual carefully before use.

After reading, please keep the manual properly for future reference.

Thank you for purchasing this product.

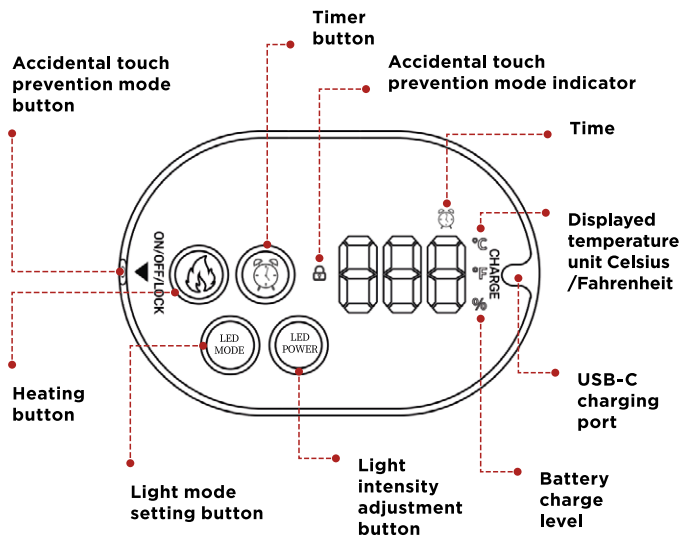
PRODUCT DESCRIPTION

This heated knee brace is light, soft and breathable, comfortable to wear and safe to use. Our knee brace is ergonomically designed with adjustable straps for knees, elbows and shoulders. It can improve blood circulation, relax muscles, relieve physical fatigue, joint pain and discomfort.

PRODUCT PARAMETERS

| | |
|--------------------------|--|
| Product Name: | BLUETOUCH - light therapy for knees |
| Working Voltage: | 5 V |
| Charging Current: | 1500 mA |
| Charging Time: | 2 hours |
| Battery Capacity: | 5000 mAh |

BUTTON DESCRIPTION

NOTE

The numbers displayed on the controller screen indicate the heating temperature ($^{\circ}\text{C}/^{\circ}\text{F}$), charging progress (%) and the set time (min). When the product is in use (not charging), the displayed number indicates the heating temperature. During use, if the anti-mistouch mode is activated, $^{\circ}\text{C}$ lights up and the number shows the current power. After a few seconds, the display returns to % and then shows the current heating temperature again. When the timer function is used, the green m indicator lights up and the number shows the set time. When the product is in charging mode, the number shows the charging progress. After being fully charged, the display will show "100 %".

OPERATION GUIDE

1. Power On: Press and hold the power button to turn on the device. If the device cannot be turned on, it may be due to insufficient battery power – please try again after charging.

2. Wear: Secure the heating brace to the knee / elbow / shoulder.

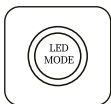
3. Heating function: Press the 🔥 button to activate the heating function. The system will automatically set to the first level: **45 °C**. There are **4 heating levels available: 45 °C - 50 °C - 55 °C - 60 °C**. Repeatedly pressing the 🔥 button switches between heating levels. The temperatures cycle from the lowest to the highest. **Note:** After the heating function is activated, the temperature increases gradually to the set value. Once the device reaches the set temperature, any further increase or decrease stops.

4. Timer function: Press the 🕒 button to activate the first timer level: **15 minutes**. The timer has **3 levels: 20 minutes - 30 minutes - 45 minutes**. Repeatedly pressing the 🕒 button switches between 30 and 45 minutes. On the fourth press, the display shows **0 min**, which means the timer function is turned off. To turn the timer on again, press the 🕒 button once more.

5. Light intensity modes: Press the LED POWER button to activate the red light function. There are four levels available: L1-L2-L3-L4. In the pulsed modes L1-L3, the light intensity gradually increases and the temperature increases gradually at the same time. L4 is a continuous mode with the highest brightness. We recommend starting with the lowest intensity and increasing it gradually.

6. Light modes: Press the LED MODE button to activate the red light mode, and "ALL" will appear on the display. Repeatedly pressing the button switches between modes: 660-850-ALL-OFF. When "OFF" is displayed, the red light mode is turned off. In ALL mode with high light intensity (L3/L4), we recommend using the device for a maximum of 30 minutes to avoid discomfort caused by excessive use. Mode descriptions: 660 - 660 nm wavelength only, 850 - 850 nm wavelength only, ALL - both wavelengths simultaneously, OFF - off.

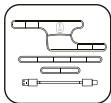
7. Anti-mistouch Touch Mode: After setting the temperature level and timing level, short press the switch button, the  icon will light up and the anti-mistouch function will be turned on. At the same time, °C lights up, the number represents the current power and jumps back to % after a few seconds, then the number again represents the current heating temperature. The controller cannot be operated in this mode to prevent accidental button touches during product use. To turn off the anti-mistouch mode, short press the power button. The  icon will turn off.

8. Switch Between °C and °F: Long press the  and  buttons at the same time for two seconds. The temperature unit will switch from °C to °F.


9. Turn Off: After use, long press the power button to turn off the knee brace. (If the timing function is set, the device will automatically shut down when the time is over, without manual operation.)

10. Storage: Keep the knee brace, extension straps and charging cable together securely to avoid losing accessories.

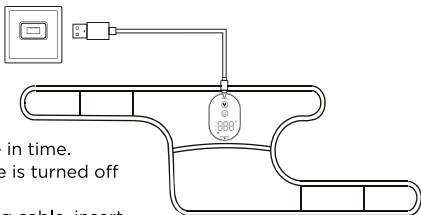CHARGING

1. When the display shows that the battery level is below 20%, please charge the device in time.
2. Make sure the knee brace is turned off before charging.
3. Use the supplied charging cable, insert it into the Type-C port of the controller, and connect the USB end to the charging device.
4. The controller display will light up, indicating that charging has started. The number shows the current battery level and will gradually increase during charging.
5. When “100%” is displayed, the knee brace is fully charged.

Note: Do not operate the controller during charging. The knee brace cannot be used while charging.

PRECAUTIONS

1. Wear thin clothing such as a shirt or yoga wear during use; thick clothing may reduce effectiveness.
2. Do not use the product directly on bare skin, as it may cause burns.
3. Prolonged use at high temperature is not recommended.
4. Do not expose the knee brace to sunlight or fire.
5. Do not disassemble or modify the heating element or control circuit.

SAFETY PRECAUTIONS

1. For first use, select a low temperature setting. Mild redness and warmth are normal. If severe redness or swelling occurs, stop using immediately and consult a doctor.
2. Avoid direct contact with the skin during long-term use.
3. Start from the lowest temperature and adjust gradually.
4. Stop using immediately if dizziness, headache, asthma, or other discomfort occurs.
5. Turn off and store the knee brace properly when not in use.
6. This product is a household appliance and not a medical device or medicine.

7. Use only according to the instructions; improper use may cause injury or property damage.

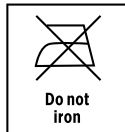
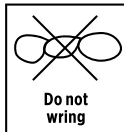
WARNING

CONSULT A DOCTOR BEFORE USE:

- Patients suffering from osteoporosis, fractures, sprains, acute pain, or other leg diseases;
- Patients with nerve disorders in the legs or sensory disturbances of the skin;
- Pregnant or postpartum women;
- Patients with heart disease;
- This product has a heating function; people who are not sensitive to temperature should pay special attention.

ABOUT CLEANING

1. Before cleaning, make sure the power is completely turned off.
2. When cleaning, wipe the knee brace with a damp cloth, being careful not to wet the controller. Do not use strong cleaners or hard brushes.
3. Lay flat to dry.
4. This product does not support machine washing or hand washing.
5. Do not dry clean to avoid damaging the heating element.

SAFETY INSTRUCTIONS

Safety Instructions

Use the device only for purposes described in the manual. Do not use if the device is visibly damaged, overheated, or emits an unusual smell. Persons with reduced heat sensitivity must use the device only under supervision. Do not use on open wounds, burned skin, fresh scars, infections, or inflamed areas. Do not use if you have an implanted electronic device (e.g., pacemaker) unless approved by a doctor. Keep out of reach of children unless stated otherwise. Do not use while sleeping. Do not expose the device to water or high humidity. Do not use near flammable materials.

Contraindications and Warnings

Not suitable for persons with photosensitive conditions. Not suitable for light-sensitive epilepsy. Not suitable during pregnancy (unless approved by a doctor). If irritation, redness, or pain occurs, stop using the device immediately.

Maintenance and Cleaning

Allow the device to cool down after use. Clean only with a soft dry or slightly damp cloth. Do not use alcohol, solvents, or aggressive cleaning agents. Do not immerse the device in water.

Storage

Store in a dry place away from direct sunlight. Storage temperature: $-10\text{ }^{\circ}\text{C}$ to $+45\text{ }^{\circ}\text{C}$. Protect from dust and humidity.

Device Disposal

Do not dispose of with household waste. Take to a specialized collection point.

Manufacturer

BlueTouch, s.r.o.

Zelený pruh 1974/109,

140 00 Prague, Czechia

ID: 09437223, VAT: CZ09437223

info@bluetouch.cz

EU COMPLIANCE

This product has been designed and manufactured in accordance with the applicable legal regulations of the European Union relating to consumer electronics.

The product complies with the requirements of the following EU directives and regulations:

- Directive 2014/30/EU on electromagnetic compatibility (EMC),
- Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS),
- ČSN EN IEC 55014-1.

Based on compliance with the above requirements, the product is marked with the CE marking.

The EU Declaration of Conformity is available from the importer and can be provided upon request to the relevant supervisory authorities or customers.

This product is intended exclusively for wellness and relaxation use.

It is not a medical device within the meaning of Regulation (EU) 2017/745 and is not intended for the diagnosis, treatment, or prevention of diseases.

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**MORE INFORMATION ABOUT LIGHT THERAPY
CAN BE FOUND ON OUR WEBSITE.**

BLUETOUCH.