REDLIGHT

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REDLIGHT

LIGHT THERAPY FOR THE MUSCULOSKELETAL SYSTEM

- **4 LIGHT MODES**
- **3 VIBRATION LEVELS**
- **3 LIGHT INTENSITY LEVELS**

WIRELESS

**UNIVERSAL
SIZE**

**5 000 mAh
BATTERY**

46 LED LIGHTS

WARNING

This device is not intended for diagnosis, treatment, or prevention of diseases. The company provides no explicit or implied guarantees regarding therapeutic effectiveness related to physical health, dermatological conditions, mental health, or any other medical field. Use of this product for therapeutic purposes must be supervised by qualified healthcare professionals and follow their instructions. Users should fully understand the intended purpose of this product. Any legal responsibility arising from misuse, health risks, or subsequent claims lies solely with the responsible party.

USER MANUAL

BLUETOOTH – Light Therapy for the Musculoskeletal System

Please read this manual carefully before use.

After reading, please keep the manual properly for future reference.

PRODUCT DESCRIPTION

Welcome and thank you for choosing our **Multifunctional Red Light Heated Massage Pad**, an innovative product that combines advanced red light technology with a comfortable massage experience. It is designed to provide physical and mental relaxation and relief through deeply penetrating red light irradiation and precise, comfortable massage. The multifunctional red light heated massage pad is suitable for various scenarios, offering a convenient massage experience at home, in the office, or while traveling. It is the ideal choice to relieve fatigue and relax both body and mind.

WARM REMINDER

This product is a household appliance and cannot replace medical equipment or medication. If any discomfort or abnormality occurs during use, please stop using it immediately and consult a professional.

PRODUCT PARAMETERS

- Material: Neoprene
- Working voltage: 5V
- Battery model: 955565
- Battery capacity: 5000mAh
- Charging current: 1850mA
- Charging time: 3.5 hours

USING STEPS

01: Before use, please check the integrity of the product to ensure there is no breakage or obvious damage. Choose a quiet and comfortable environment so you can fully relax and enjoy the massage.

02: Wear the product and extension strap on the area to be used, and adjust the position to ensure it fits tightly on the massage area.

03: Long press the power button on the controller to turn on the device. If it cannot be turned on normally, it may be due to insufficient battery power. Please charge it and try again.

04: Touch the vibration button to turn on the massage function. The display shows “1”, indicating vibration level 1. Touch the vibration button repeatedly to switch vibration levels: 1–2–3–0. The display shows “0” when the massage function is turned off. *When using for the first time, it is recommended to start with a lower level and gradually increase the intensity.*

05: Touch the red light mode button to turn on the red light mode, and “660” will be displayed. Touch the button repeatedly to switch modes: 660–850–F10–ALL–OFF. The display shows “OFF” when the red light mode is turned off. *It is recommended not to use for more than 30 minutes per session.* 660 – 660 nm wavelength only, 850 – 850 nm wavelength only, F10 – alternating wavelengths 660 and 850, ALL – both wavelengths simultaneously, OFF – switched off.

06: After turning on the red light mode, touch the brightness button to adjust the red light brightness. The brightness cycles through L1–L2–L3. *If the red light mode is not turned on, skip this step. If there is no response when touching the brightness button, please check whether the red light mode is activated.*


07: Touch the timing button to turn on the timing function. The default timing is 15 minutes. Touch repeatedly to switch between 30–45–15 minutes. If no time is set, the device will automatically turn off after 15 minutes.

08: After completing the settings, short press the power button to activate the button lock. The display shows the current power, and the controller buttons are locked. To unlock, press the power button again.

09: During use, long press the power button until the display shows “OFF” to turn off the device. The device will also automatically shut down when the set time ends.

10: After use, store the product and accessories in a dry and ventilated place, away from direct sunlight and moisture.

Note: If you feel uncomfortable or notice any abnormalities during use, please stop using the product immediately.

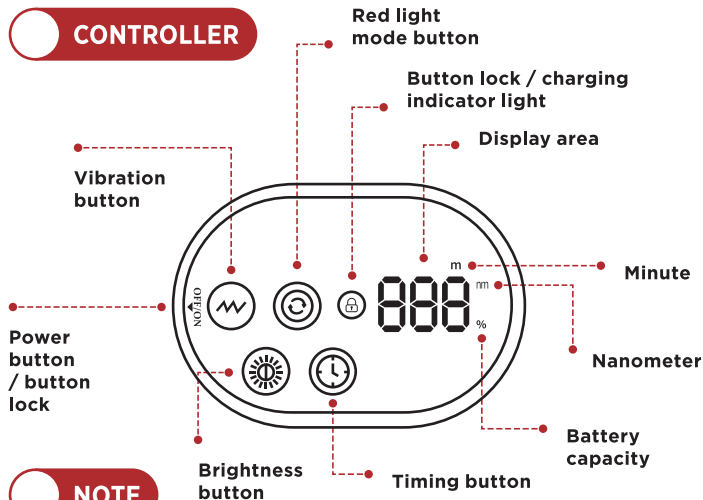The numbers displayed on the controller screen indicate the red light wavelength (nm), power and charging progress (%) and time (min).

1. When the red light mode is set, the "nm" symbol indicates the wavelength of the red light.
2. When the vibration level is set, the number indicates the massage intensity.
3. When using the timing function, the "m" symbol indicates the set time.
4. When the product is in use (not charging), the number + % indicates the current power.
5. When the product is charging, the number + % indicates the charging progress.
6. The brightness button is used to adjust the light intensity. Three intensity levels are available: L1/L2/L3. We recommend starting with the lowest intensity and gradually increasing it according to your feelings and the body's reactions.

ABOUT CHARGING

1. When the device is turned on, briefly press the power button and the display area will show the current battery level.

2. When the screen shows that the power is lower than 20%, charge the device in time.
3. Before charging, make sure that the device is turned off.
4. Take out the supplied charging cable and connect it to the controller and the charging adapter. The controller display lights up, indicating that the device has started charging, and the numbers represent the current power. During charging, the number will slowly increase.
5. The controller display shows "100%", which means that the device is fully charged.

Note: When charging, please use our original charging cable, otherwise problems such as failure to charge or incomplete charging may occur. Do not use the controller during charging, otherwise it may be damaged.

WARNINGS

- This product is a household appliance and cannot replace medical equipment or medicines.
 - When not in use, turn off the product and store it in a safe place.
 - Improper use may lead to safety accidents or property damage.
1. Before use, make sure that the product is completely dry.
 2. This product is not waterproof and must not be used in a humid environment.
 3. It is recommended that the product be in contact with the skin at zero distance during use. Wearing clothing may affect the effect of use.
 4. During operation, do not look directly at the beads.
 5. The recommended optimal single exposure time of red light is 15-30 minutes.
 6. If you feel unwell during use, stop using the device immediately.
 7. During use, avoid excessive bending or pulling of the product to prevent damage to the internal circuit.
 8. Do not disassemble or modify the product to avoid artificial damage.
 9. If a product failure or abnormal situation occurs, contact customer service in a timely manner.
 10. Do not expose the product to sunlight or fire to avoid damage to the thermal conductor of the insulator.

WARNING



CONSULT A DOCTOR BEFORE USE:

- Pregnant and breastfeeding women
- People with acute illnesses or inflammatory diseases
- Patients with heart diseases and people with pacemakers
- Patients recovering from surgery
- People allergic to red light or massage
- People with a tendency to bleed or who are bleeding
- People with malignant tumors
- People with severe neurological disorders
- Persons lacking legal capacity

Before use, it is recommended to assess the application based on your personal health condition and your doctor's advice. If you have any questions or feel uncomfortable, stop using the device immediately and consult a professional doctor.

CLEANING AND MAINTENANCE

1. Before cleaning, please make sure to disconnect the power completely.
2. When cleaning, wipe the product with a damp cloth, but be careful not to wet the Type-C interface.
3. Do not use strong cleaners or hard brushes.
4. Lay flat to dry.
5. This product does not support machine wash or hand wash.
6. Do not dry clean.

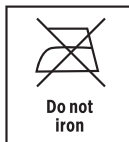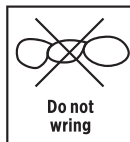WARRANTY & AFTER-SALES SERVICE

Thank you for purchasing our product. Your satisfaction and support are our greatest motivation. Please contact us directly if you have any questions. We will provide prompt service if you encounter any problem with our product. **Thank you!**

SAFETY INSTRUCTIONS

Safety Instructions

Use the device only for purposes described in the manual. Do not use if the device is visibly damaged, overheated, or emits an unusual smell. Persons with reduced heat sensitivity must use the device only under supervision. Do not use on open wounds, burned skin, fresh scars, infections, or inflamed areas. Do not use if you have an implanted electronic device (e.g., pacemaker) unless approved by a doctor. Keep out of reach of children unless stated otherwise. Do not use while sleeping. Do not expose the device to water or high humidity. Do not use near flammable materials.

Contraindications and Warnings

Not suitable for persons with photosensitive conditions. Not suitable for light-sensitive epilepsy. Not suitable during pregnancy (unless approved by a doctor). If irritation, redness, or pain occurs, stop using the device immediately.

Maintenance and Cleaning

Allow the device to cool down after use. Clean only with a soft dry or slightly damp cloth. Do not use alcohol, solvents, or aggressive cleaning agents. Do not immerse the device in water.

Storage

Store in a dry place away from direct sunlight. Storage temperature: -10 °C to +45 °C. Protect from dust and humidity.

Device Disposal

Do not dispose of with household waste. Take to a specialized collection point.

Manufacturer

BlueTouch, s.r.o.
Zelený pruh 1974/109,
140 00 Prague, Czechia
ID: 09437223, VAT: CZ09437223
info@bluetouch.cz

EU COMPLIANCE

This product has been designed and manufactured in accordance with the applicable legal regulations of the European Union relating to consumer electronics.

The product complies with the requirements of the following EU directives and regulations:

- Directive 2014/30/EU on electromagnetic compatibility (EMC),
- Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS),
- ČSN EN IEC 55014-1.

Based on compliance with the above requirements, the product is marked with the CE marking.

The EU Declaration of Conformity is available from the importer and can be provided upon request to the relevant supervisory authorities or customers.

This product is intended exclusively for wellness and relaxation use.

It is not a medical device within the meaning of Regulation (EU) 2017/745 and is not intended for the diagnosis, treatment, or prevention of diseases.

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**MORE INFORMATION ABOUT LIGHT THERAPY
CAN BE FOUND ON OUR WEBSITE.**

BLUETOUCH.