REDLIGHT

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BEZDRÁTOVÉ

KOMPAKTNÍ

**BATERIE
1 000 mAh**

3 LED DIODY

BLUETOOTH.

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1. PRODUCT DESCRIPTION & USER MANUAL

BLUETOUCH – Light Therapy; Skin Care

Thank you for choosing BLUETOUCH – Light Therapy; Skin Care. This beauty device combines modern skincare technologies, including LED light therapy, gentle vibrations, heating, ion technology, and EMS microcurrent. It is designed for your at-home beauty routine to help promote a fresher, smoother, and more radiant-looking complexion. The device offers the following four treatment modes.

2. MODE DESCRIPTION

1. REPAIR – A mode designed for gentle skin regeneration. It combines vibrations, heating, red LED light, and infrared light. It helps relax and gently warm the skin, leaving it looking refreshed and well-rested. It is particularly suitable as a relaxing step in your evening skincare routine or whenever your skin appears tired. Technologies used: vibrations, heating, red LED light, infrared light.

2. CLEAN – A mode designed to promote a cleaner and fresher-looking complexion. It combines vibrations, heating, negative ions, blue LED light, and red LED light. We recommend using it after removing makeup and cleansing the skin. It helps prepare the skin for the next steps of your beauty routine. Technologies used: vibrations, heating, negative ions, blue LED light, red LED light.


3. ABSORB – Designed for use with facial serum, cream, or other skincare products. It combines vibrations, heating, EMS microcurrent, blue LED light, red LED light, and infrared light. It helps evenly distribute skincare products into the skin while enhancing the feeling of intensive skincare. After use, the skin may appear smoother, fresher, and better nourished. Technologies used: vibrations, heating, EMS microcurrent, blue LED light, red LED light, infrared light.

4. LIFTING – This mode is designed to promote a firmer and more lifted-looking complexion. It combines vibrations, heating, positive ions, blue LED light, and infrared light. It is suitable for regular care of the facial contours, cheeks, chin, and other areas where the skin may appear loose or tired. Technologies used: vibrations, heating, positive ions, blue LED light, infrared light.


3. PRODUCT SPECIFICATIONS

Handheld beauty device

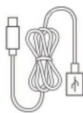USB cable

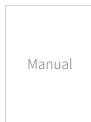User manual

Product Name	BLUETOUCH – Light Therapy; Skin Care
Product Dimensions	82.6 × 47 × 50 mm
Output Power	2 W
Light Wavelengths	850 nm, 630 nm, 465 nm
Treatment Modes	repair, clean, absorb, lifting
Battery Capacity	1000 mAh
Charging Port	USB-C
Charging Input	DC 5 V / 1 A

4. FUNCTIONS & OPERATION

Powering the Device On and Off – To switch on the device, press and hold the power button for 1 second. Briefly press the function button to switch between the individual treatment modes. The device offers four modes. If no mode is selected after switching on, the device will automatically start in Repair mode at Intensity Level 1. The device automatically switches off after 5 minutes of operation. An audible signal will sound before the device powers off.

Recommendations for Use – Apply the device gently to the skin without excessive pressure. Move it slowly in upward strokes, working from the centre of the face outwards. When using Absorb mode, apply a facial serum or cream to the skin and then use the device to distribute the product evenly. Do not keep the device in one place for too long. Each treatment mode is set to run for 5 minutes.

Mode 1: Repair – After switching on the device and placing the metal treatment head in contact with the skin, the device automatically activates Repair mode. The corresponding mode indicator and intensity level indicator will illuminate.

Mode 2: Clean – When the device is in Repair mode, briefly press the power button. The device will emit one beep and switch to Clean mode. The blue operating indicator and the intensity level indicator will illuminate.

Mode 3: Absorb – Using the same sequence, the device switches from Clean mode to Absorb mode. The corresponding operating indicator and intensity level indicator will illuminate.

Mode 4: Lifting – When the device is in Absorb mode, briefly press the power button. The device will emit one beep and switch to Lifting mode. The operating indicator and intensity level indicator will illuminate.

5. INTENSITY SETTINGS

The + and – buttons are used to adjust the vibration intensity.

– Button

Briefly press the – button to decrease the vibration intensity. A beep will sound, and the device will switch to a lower intensity level.

+ Button

Briefly press the + button to increase the vibration intensity. After the first press, the device switches to Intensity Level 2. Two white indicator lights will illuminate.

After the second press, the device switches to Intensity Level 3. Three white indicator lights will illuminate.

After the third press, the device switches to Intensity Level 4. Four white indicator lights will illuminate.

6. RECOMMENDED USE

Before use, thoroughly remove makeup and cleanse the skin. Select the treatment mode that best suits your skin's current needs and gently glide the device over your face. When using Absorb mode, apply a facial serum or cream and move the device slowly in upward strokes, working from the centre of the face outwards. The recommended treatment time is a few minutes per beauty routine. Do not keep the device in one place for too long. Guide the device gently without applying excessive pressure:

- from the chin towards the jawline,
- from the corners of the mouth towards the ears,
- from the nose across the cheeks outwards,
- from the centre of the forehead towards the temples,
- along the neck in upward strokes.

During use, the skin should feel comfortably firm, but it should not become irritated.

Do not use the device on injured, irritated, or inflamed skin.

Do not use it around the eyes, on the eyelids, or directly on mucous membranes. This device is not a medical device.

It is intended solely for cosmetic skincare purposes. Individuals with a pacemaker, metal facial implants, epilepsy, serious skin conditions, those who are pregnant, or anyone with medical concerns should consult a physician before use. If you experience discomfort, burning, pain, or significant skin irritation, stop using the device immediately.

7. DISCLAIMER

1. This product is not intended to diagnose, treat, cure, or prevent any disease. The company makes no claims, representations, or warranties regarding the ability of this product to treat any physical, skin, or mental condition. If you have any health condition requiring medical attention, consult a qualified healthcare professional.
2. The company accepts no legal liability for any claims arising from safety issues caused by improper use of the product.

8. EU REGULATORY COMPLIANCE

This product has been designed and manufactured in accordance with the applicable European Union legislation governing consumer electronic products.

The product complies with the requirements of the following EU directives and regulations:

- Directive 2014/30/EU on Electromagnetic Compatibility (EMC),

- Directive 2011/65/EU on the Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment (RoHS),
- ČSN ISO 13485
- ČSN EN IEC 55014-1

Based on compliance with the above requirements, this product bears the CE marking.

The EU Declaration of Conformity is available from the importer and can be provided upon request to the relevant supervisory authorities or customers.

This product is intended exclusively for wellness and relaxation purposes. It is not a medical device within the meaning of Regulation (EU) 2017/745 and is not intended to diagnose, treat, cure, or prevent any disease.

9. SAFETY INSTRUCTIONS

Safety Instructions

Use the device only for the purposes described in this manual.

Do not use the device if it is visibly damaged, overheated, or emits an unusual odour. Persons with reduced sensitivity to heat should use the device only under supervision.

Do not use on open wounds, burned skin, fresh scars, infections, or inflamed areas. Do not use the device if you have an implanted electronic device (e.g. a pacemaker), unless otherwise advised by your physician. Keep out of reach of children unless otherwise stated. Do not use the device while sleeping. Do not expose the device to water or excessive humidity. Do not use the device near flammable materials.

Contraindications and Warnings

Not suitable for individuals with photosensitive conditions.
Not suitable for individuals with photosensitive epilepsy.
Not suitable during pregnancy unless approved by a physician.
If irritation, redness, or pain occurs, discontinue use immediately.

Maintenance and Cleaning

Allow the device to cool down after use.
Clean only with a soft, dry, or slightly damp cloth.
Do not use alcohol, solvents, or aggressive cleaning agents.
Do not immerse the device in water.

Storage

Store in a dry place away from direct sunlight.
Storage temperature: -10 °C to +45 °C.
Protect from dust and moisture.

Disposal

Do not dispose of the device with household waste.
Take it to a designated electronic waste collection point.

Manufacturer

BlueTouch, s.r.o.
Zelený pruh 1974/109,
140 00 Prague, Czechia
ID: 09437223, VAT: CZ09437223
info@bluetouch.cz