REDLIGHT

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BEZDRÁTOVÉ

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**BATERIE
2 000 mAh**

280 LED DIOD

OBSAH

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 **BLUETOUCH.**

CONTENTS

1. Specification details	55
2. Mode switching.	56
3. Charging the mask.	56
4. Light therapy modes	57
5. Recommended frequency of use	57
6. How to use.	58
7. Step by step	59
8. Package contents.	60
9. Safety instructions.	61
10. Danger	61
11. Caution.	62
12. Caution.	63
13. Cleaning and maintenance.	63
14. Disclaimer	64
15. EU compliance	64
16. Safety instructions.	65

1. SPECIFICATION DETAILS


Product name	Redlight light therapy for the face
Model	For the face
Product dimensions	634 × 228 × 154 mm
Product weight	318 g
Output power	4.2 W
Number of LEDs	280 LED diodes
9 light modes	R+N, red, NIR, blue, green, purple, turquoise, yellow, white
Near infrared light	850 nm
Red light	630 nm
Blue light	465 nm
Green light	530 nm
Purple light	blue light 465 nm + red light 630 nm
Indigo light	blue light 465 nm + green light 530 nm
Yellow light	green light 530 nm + red light 630 nm
White light	blue light 465 nm + green light 530 nm + red light 630 nm
Battery capacity	2000 mAh
Charging port	USB-C
Charging input	DC 5 V / 1 A
Charging time	3.5 h
Operating time	1.5 h

2. MODE SWITCHING

Press the **M button** on the remote control to switch between modes:

Remote control

- M1** R+N (red + infrared)
- M2** RED - red light
- M3** NIR - near infrared light
- M4** BLU - blue light
- M5** GRN - green light
- M6** PURP - purple light
- M7** TEAL - turquoise light
- M8** YEL - yellow light
- M9** WHT - white light

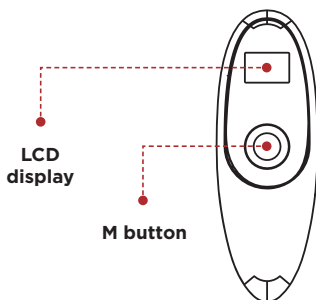3. CHARGING THE MASK

Connect the device to a power source using a **USB-C cable**.

Charging indicators

Red light: charging in progress

Green light: fully charged


Important notes

Charging time: approx. 3.5 hours

Battery life: up to 90 minutes of use

The mask **cannot be used while charging**.

When the battery level is low, the indicator light on the upper part of the mask will turn on. It is recommended to recharge the device in time so that further use is not affected.

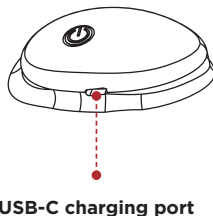4. LIGHT THERAPY MODES

M1 Red + near infrared light

Benefits: helps reduce wrinkles and fine lines, supports skin microcirculation.

M2 Red light

Benefits: stimulates cell activity and supports skin metabolism.

M3 Near infrared light (NIR)

Benefits: helps regenerate damaged skin, supports firmness and elasticity.

M4 Blue light

Benefits: helps eliminate bacteria and supports skin blood circulation.

M5 Green light

Benefits: soothes the skin and helps unify its tone.

M6 Purple light

Benefits: supports the production of collagen and elastin, contributes to firmer skin.

M7 Turquoise light

Benefits: increases the energy of skin cells and accelerates metabolism.

M8 Yellow light

Benefits: supports skin regeneration, especially suitable after fatigue or sun irritation, helps brighten the skin.

M9 White light

Benefits: penetrates into deeper layers of the skin and increases the overall effect of the treatment.

5. RECOMMENDED FREQUENCY OF USE

- use **3-4× per week 10-20 minutes** per one procedure
- with good tolerance, the care can be gradually intensified as needed

6. HOW TO USE

There are two ways to control the mask:

1. Control by the button on the mask

Press and hold the power button for 2 seconds – a short sound signal will be heard and the mask will enter standby mode, while the green indicator starts flashing.

By briefly pressing the button, you switch the light modes **M1 to M9**. Each mode has a **preset time of 10 minutes**, after which the mask switches off automatically.


By pressing and holding the button again for 2 seconds you switch the mask off.

2. Control using the remote control

Using the remote control, it is possible to switch the light modes **M1 to M9**. The remote control **cannot be used to turn the mask on or off**.

Tips for safe use

- Before first use, test the mask **on your arm for 1 minute**.
- If no unpleasant sensations appear, you can use the product normally.
- If the device needs to be restarted, simply turn it on again.
- If you are not using the mask, **switch it off and store it back in the protective packaging**.

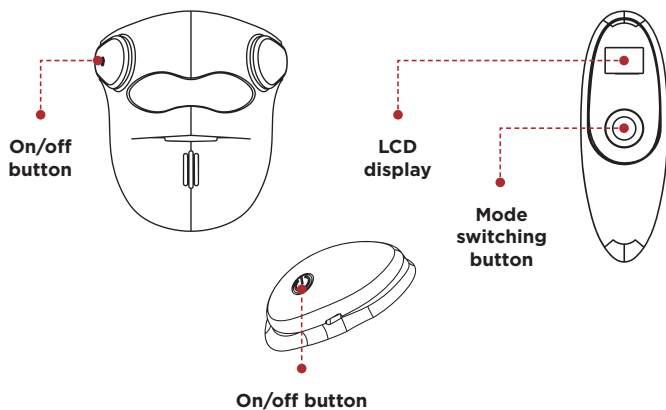
7. STEP BY STEP

- Step 1:** Clean the skin, the skin should be clean and dry.
- Step 2:** Adjust the strap length and put on the face mask so that it fits comfortably.
- Step 3:** Turn on the device by pressing and holding the power button – the light will turn on.
- Step 4:** Briefly press the **M button** controller to switch between modes **M1 to M9**.
- Step 5:** One procedure lasts **10 minutes**.
- Step 6:** Enjoy the light therapy – after the program ends, the device **switches off automatically**.

After the treatment

After completing the procedure, you can **hydrate the skin according to your usual skincare routine** for maximum effect.


8. PACKAGE CONTENTS

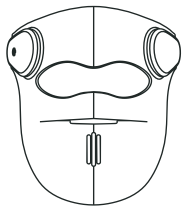
LED mask


Remote control

User manual

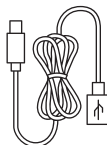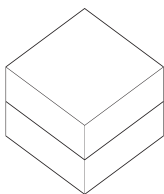USB-C charging cable

Box

9. SAFETY INSTRUCTIONS

For safe use of this device, carefully read the manual and follow the instructions given.

-  **Danger:** situations that may lead to serious injury or danger to life.
-  **Warning:** situations that may cause serious injury to persons.
-  **Caution:** situations that may cause injury or material damage.

10. DANGER

Persons using medical devices such as **pacemakers, ECG devices, heart and lung support devices, or infusion pumps** should not use this product. Malfunction of the medical device or physical discomfort could occur.

Do not use the product near a **filled bathtub, shower, sink, or other containers with water**. Keep the device and the power supply dry.

Do not use the device near a **stove, heating, or air conditioning units**, to avoid damage to the product or an accident.

If the **product, remote control, or power cable is damaged**, do not use it – there is a risk of electric shock.

Keep the device out of reach of **infants and children**, so that dangerous situations do not occur.

The product **is not waterproof**. Do not immerse it in water, as damage or an accident could occur. If it gets wet, do not continue using the device.

11. CAUTION

- This product is intended **only for use on the face**. Do not use it to solve problems other than those for which it is intended.
- **Do not disassemble, modify, or repair the device yourself** to ensure its proper function.
- Do not use simultaneously **with photosensitive (light-sensitive) ingredients or cosmetics containing acids (e.g. AHA, BHA, retinoids, etc.)**.

The product is not suitable for the following persons:

- persons **under 18 years of age**
- patients with **heart disease**
- persons with **abnormal blood pressure**
- persons with **mental illness**
- persons with **increased sensitivity to light (photosensitivity)**
- persons with **an infectious disease**
- persons with **diabetes**
- persons with **melanoma or other skin tumors**
- persons with **malignant tumors**
- **pregnant and breastfeeding women** and persons with very sensitive skin
- skin **damaged by the sun**
- persons taking medication, especially **oncology patients** or persons with **precancerous moles**
- persons with **wound healing disorders**
- persons with **diseases of the immune system**
- persons who have undergone **botox application in the last 3 months**
- persons with **hyperpigmentation or contact dermatitis** caused by chemical substances
- persons taking **steroid hormonal medication** long-term or suffering from **telangiectasia** (dilated vessels) caused by a liver disorder
- persons with **dermal fillers** – consult a doctor before use
- if you are currently undergoing **medical treatment**, consult use with a doctor in advance

Areas unsuitable for use

Using the product on the following areas may cause undesirable reactions or skin problems:

- areas other than the intended ones
- areas irritated by cosmetics
- inflamed areas, for example acne
- skin damaged by sunlight
- age spots (solar lentigo)
- more extensive injuries or open wounds
- areas with irritation caused by itching, redness, or mechanical irritation

12. CAUTION

- Keep your **hands dry** when using the device.
- The product is a **precision electronic device** – do not bend it, pull it, or sit on it to avoid damage or short circuit. Store it properly after use.
- Follow the recommended **time and frequency of use** to avoid unpleasant feelings.
- If you feel **any discomfort during use, stop using the device immediately.**
- Protect both the main unit and the power supply from **impacts and damage.**

13. CLEANING AND MAINTENANCE

To maintain maximum effectiveness and long service life of the device: After each use, store the product in the protective box so that it is not damaged or contaminated with dust. If necessary, clean the device with a **soft lint-free cloth** or non-woven fabric slightly moistened with alcohol – wipe the parts where dirt is present.

Product disposal

This product contains a battery. **Do not throw it into regular household waste.** Hand it over at the appropriate collection point for recycling of electronic waste.

14. DISCLAIMER

1. This product **is not intended** for the diagnosis, treatment, or prevention of any disease. The company makes no claims, statements, or warranties regarding the ability of this product to treat any physical, dermatological, or mental condition. If you have any health issues requiring medical attention, please consult a qualified healthcare professional.
2. The company assumes no legal liability for any claims arising from safety issues caused by **improper** use of the product.

15. EU COMPLIANCE

This product has been designed and manufactured in accordance with the applicable legal regulations of the European Union relating to consumer electronics.

The product complies with the requirements of the following EU directives and regulations:

- Directive 2014/30/EU on electromagnetic compatibility (EMC),
- Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS),
- ČSN ISO 13485,
- ČSN EN IEC 55014-1.

Based on compliance with the above requirements, the product is marked with the CE marking.

The EU Declaration of Conformity is available from the importer and can be provided upon request to the relevant supervisory authorities or customers.

This product is intended exclusively for **wellness and relaxation use**.

It is not a medical device within the meaning of Regulation (EU) 2017/745 and is not intended for the diagnosis, treatment, or prevention of diseases.

16. SAFETY INSTRUCTIONS

Safety Instructions

Use the device only for purposes described in the manual. Do not use if the device is visibly damaged, overheated, or emits an unusual smell. Persons with reduced heat sensitivity must use the device only under supervision. Do not use on open wounds, burned skin, fresh scars, infections, or inflamed areas. Do not use if you have an implanted electronic device (e.g., pacemaker) unless approved by a doctor. Keep out of reach of children unless stated otherwise. Do not use while sleeping. Do not expose the device to water or high humidity. Do not use near flammable materials.

Contraindications and Warnings

Not suitable for persons with photosensitive conditions. Not suitable for light-sensitive epilepsy. Not suitable during pregnancy (unless approved by a doctor). If irritation, redness, or pain occurs, stop using the device immediately.

Maintenance and Cleaning

Allow the device to cool down after use. Clean only with a soft dry or slightly damp cloth. Do not use alcohol, solvents, or aggressive cleaning agents. Do not immerse the device in water.

Storage

Store in a dry place away from direct sunlight.
Storage temperature: -10 °C to +45 °C. Protect from dust and humidity.

Device Disposal

Do not dispose of with household waste.
Take to a specialized collection point.

Manufacturer

BlueTouch, s.r.o.
Zelený pruh 1974/109,
140 00 Prague, Czechia
ID: 09437223, VAT: CZ09437223
info@bluetouch.cz