REDLIGHT

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REDLIGHT

LIGHT THERAPY – HANDHELD LAMP

- 2 LIGHTING MODES
- 4 PRESET MODES
- 5 WAVELENGTHS

WIRELESS

COMPACT

BATTERY
1 800 mAh

5 LEDS

CONTENTS

1. Product Specifications	36
2. Function Instructions	37
3. Digital Display Functions	38
4. Device Maintenance	38
5. Disclaimer	39
6. EU Compliance	39
7. Safety Instructions	40

1) Although many people experience pain relief almost immediately after starting red light therapy, for others measurable relief develops gradually over time. The body's response to therapy varies depending on the type and severity of the condition. Therefore, we recommend using red light therapy regularly and over the long term. If you notice results immediately, that's great—continue using it as needed. If you do not experience relief right away, please be aware that this is also completely normal, so do not be discouraged. Simply treat the affected area regularly with three short treatments per day, each lasting up to 5 minutes per treatment area. After one month of consistent red light therapy, most people notice positive effects.

2) Before starting red light therapy, consult your physician and do not discontinue any other prescribed treatments. Allow the device sufficient time to provide optimal care. Hold the device directly against the affected area or a few centimetres above it, depending on the desired treatment depth and coverage. Red light penetrates several centimetres into tissues and joints.

DO NOT OPERATE THE DEVICE FOR MORE THAN 5 MINUTES AT A TIME.

If you wish to treat multiple areas of the body consecutively, switch the device off for at least 1 minute after every 5 minutes of operation. This helps prevent the unit from overheating.

Treat each treatment area only once during a single treatment session.


1. PRODUCT SPECIFICATIONS

Parameter	Value	Parameter	Value
Product Dimensions	154 × 34,2 mm (without attachment)	Net Weight	0,17 kg
Material	Stainless Steel	Wavelength	460 / 630 / 660 / 850 / 900 nm
Irradiance	Surface: 217 mW/cm ²	Irradiance	at a Distance of 7.5 cm: 51 mW/cm ²
Illuminance	180 lux at 15 cm	Battery	Li-Ion 18650


630 + 660 nm Red Light

Relieves pain, promotes
wound healing

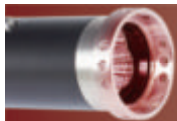850 + 900 nm Near-Infrared Light


Deep tissue penetration for enhanced
treatment

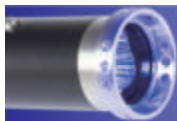460 nm Blue Light

Helps cleanse the skin
and eliminate acne-causing bacteria

2. FUNCTION INSTRUCTIONS

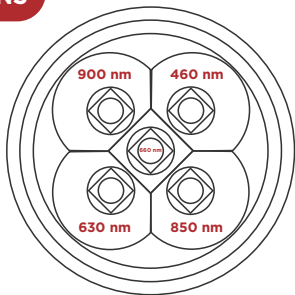1. Timer:

The default setting is 5 minutes. Alternatively, the timer can be disabled by pressing and holding the button for 5 seconds.

2. Four Modes:

- (a) **Mode 1 – All LEDs active:**
460 nm, 630 nm, 660 nm, 850 nm, 900 nm
- (b) **Mode 2 – Red light only:**
630 nm, 660 nm
- (c) **Mode 3 – Near-infrared light only:**
850 nm, 900 nm (barely visible to the naked eye)
- (d) **Mode 4 – Blue light only:**
460 nm

3. Supports pulse mode.

4. Button Operation:

1st press: Mode 1

2nd press: Mode 1 + 10 Hz pulse mode

3rd press: Mode 2

4th press: Mode 2 + 10 Hz pulse mode

5th press: Mode 3 (barely visible to the naked eye)

6th press: Mode 3 + 10 Hz pulse mode

7th press: Mode 4

8th press: Mode 4 + 10 Hz pulse mode

9th press: Power off

10th press: The lamp can be switched off at any time by pressing the button. When switched on again, the device remembers the last active mode (memory function).

Research indicates a brief “resting phase” that occurs within cells during pulsed near-infrared (NIR) light therapy. Although it lasts only a few milliseconds, this short pause may have a significant effect by allowing cells under oxidative stress to respond more effectively to light. As a result, pulsed therapy may provide greater stimulation of cellular activity and support tissue regeneration.

3. DIGITAL DISPLAY FUNCTIONS

1. Battery level indicator
2. 5-minute countdown with automatic shut-off after 5 minutes.


M1:
All LEDs active
+ Flashing

M2:
Red light only
+ Flashing

M3:
Near-infrared
light only
+ Flashing

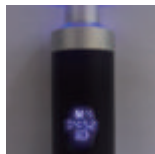M4:
Blue light only
+ Flashing

Warning

Do not use the device near water. Do not look directly into the light. Always wear protective eyewear when using the device. Do not perform treatment near the eyes.

4. DEVICE MAINTENANCE

Cleaning and Disinfection

Use only a dry cloth to clean the exterior of the device. Never use abrasive cleaning agents or immerse the device in water or any other liquid. This device must not be sterilised using liquids or an autoclave.

Storage

Store the device in a cool, dry place away from direct sunlight and excessive heat.

Warning

Prevent moisture from coming into contact with the housing or the light-emitting part of the device. Do not immerse the device in water.

Warning

If the device will not be used for an extended period (1 month or longer), remove the battery to prevent corrosion.

5. DISCLAIMER

1. This product **is not intended** for the diagnosis, treatment, or prevention of any disease. The company makes no claims, statements, or warranties regarding the ability of this product to treat any physical, dermatological, or mental condition. If you have any health issues requiring medical attention, please consult a qualified healthcare professional.
2. The company assumes no legal liability for any claims arising from safety issues caused by **improper** use of the product.
3. **Built-in Rechargeable Battery**
This product contains a built-in rechargeable battery. Do not disassemble, puncture, or expose the battery to fire, water, high temperatures, or mechanical damage. Do not charge the product unattended, on a flammable surface, or near flammable materials. If the product overheats, emits an odor, changes shape, or shows signs of damage, stop using it immediately and contact customer service.

6. EU COMPLIANCE

This product has been designed and manufactured in accordance with the applicable legal regulations of the European Union relating to consumer electronics.

The product complies with the requirements of the following EU directives and regulations:

- Directive 2014/30/EU on electromagnetic compatibility (EMC),
- Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS),
- ČSN EN IEC 55014-1.

Based on compliance with the above requirements, the product is marked with the CE marking. The EU Declaration of Conformity is available from the importer and can be provided upon request to the relevant supervisory authorities or customers. This product is intended exclusively for **wellness and relaxation use**. It is not a medical device within the meaning of Regulation (EU) 2017/745 and is not intended for the diagnosis, treatment, or prevention of diseases.

7. SAFETY INSTRUCTIONS

Safety Instructions

Use the device only for purposes described in the manual. Do not use if the device is visibly damaged, overheated, or emits an unusual smell. Persons with reduced heat sensitivity must use the device only under supervision. Do not use on open wounds, burned skin, fresh scars, infections, or inflamed areas. Do not use if you have an implanted electronic device (e.g., pacemaker) unless approved by a doctor. Keep out of reach of children unless stated otherwise. Do not use while sleeping. Do not expose the device to water or high humidity. Do not use near flammable materials.

Contraindications and Warnings

Not suitable for persons with photosensitive conditions. Not suitable for light-sensitive epilepsy. Not suitable during pregnancy (unless approved by a doctor). If irritation, redness, or pain occurs, stop using the device immediately.

Maintenance and Cleaning

Allow the device to cool down after use. Clean only with a soft dry or slightly damp cloth. Do not use alcohol, solvents, or aggressive cleaning agents. Do not immerse the device in water.

Storage

Store in a dry place away from direct sunlight.
Storage temperature: -10 °C to +45 °C. Protect from dust and humidity.

Device Disposal

Do not dispose of with household waste.
Take to a specialized collection point.

Manufacturer

BlueTouch, s.r.o.

**Zelený pruh 1974/109,
140 00 Prague, Czechia**

ID: 09437223, VAT: CZ09437223

info@bluetouch.cz