# REDLIGHT

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100 LED DIODES

### OPERATING INSTRUCTIONS

**1. Power on:** Press and hold the main button for three seconds to turn on the device. If the device does not turn on normally, charge it before use. Do not use the device while charging.

**2. Temperature setting:** Click the temperature button in the lower left corner to switch temperature levels. Available settings: 45°C / 50°C / 55°C.

**3. Red light mode:** Click the middle button to switch red-light intensity. There are 5 levels: L1 / L2 / L3 / L4 / L5. The light brightness gradually increases at levels L1-L3, and the temperature increases at the same time. L4 / L5 are pulsed modes. We recommend starting with the lowest intensity and gradually increasing it.

**4. Timer setting:** Click the timer button in the lower right corner to switch between 10 / 20 / 30 minutes.

**5. Lock (accidental-touch protection):** Press the main button to activate the lock mode. A lock icon will appear on the display. Press again to unlock.

**6. Power off:** When not in use, press and hold the main button for three seconds to turn off the device.


## BASIC INFORMATION

|                                    |                                                    |
|------------------------------------|----------------------------------------------------|
| <b>Product name:</b>               | <b>BLUETOUCH - light therapy for neck and nape</b> |
| <b>LED diodes:</b>                 | <b>100 pcs</b>                                     |
| <b>Voltage / Battery capacity:</b> | <b>5V 2A / 5 000 mAh</b>                           |
| <b>Wavelengths:</b>                | <b>660 nm + 850 nm</b>                             |
| <b>Brightness levels:</b>          | <b>L1 / L2 / L3 / L4 / L5</b>                      |
| <b>Timer:</b>                      | <b>10 / 20 / 30 minutes</b>                        |
| <b>Temperature:</b>                | <b>45°C / 50°C / 55°C</b>                          |

## WARNING

This device is not intended for diagnosis, treatment, or prevention of diseases. The company provides no explicit or implied guarantees regarding therapeutic effectiveness related to physical health, dermatological conditions, mental health, or any other medical field. Use of this product for therapeutic purposes must be supervised by qualified healthcare professionals and follow their instructions. Users should fully understand the intended purpose of this product. Any legal responsibility arising from misuse, health risks, or subsequent claims lies solely with the responsible party.

## ADDITIONAL INFORMATION
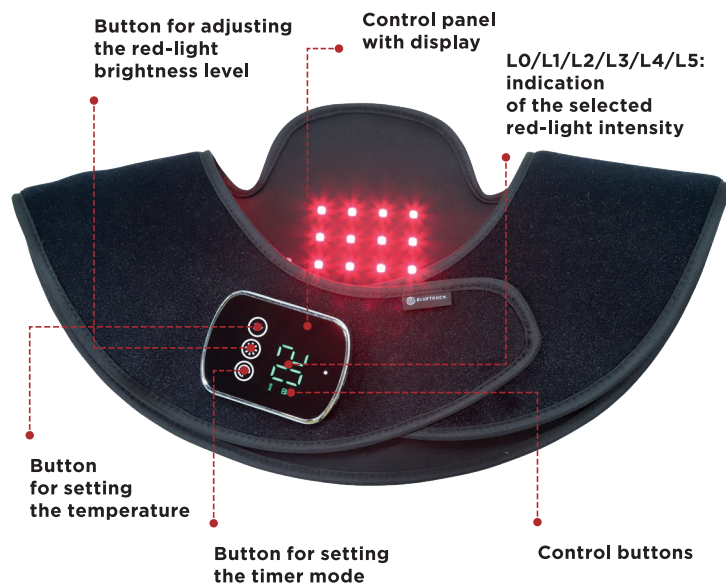
## ADDITIONAL INFORMATION

- The device automatically switches off after 30 minutes of use. To use it again, turn it back on.
- Do not use the device while it is charging.
- The device charges via a USB-C port. A charging cable is included.

## SAFETY INSTRUCTIONS

### Safety Instructions

Use the device only for purposes described in the manual. Do not use if the device is visibly damaged, overheated, or emits an unusual smell. Persons with reduced heat sensitivity must use the device only under supervision. Do not use on open wounds, burned skin, fresh scars, infections, or inflamed areas. Do not use if you have an implanted electronic device (e.g., pacemaker) unless approved by a doctor. Keep out of reach of children unless stated otherwise. Do not use while sleeping. Do not expose the device to water or high humidity. Do not use near flammable materials.

### Contraindications and Warnings

Not suitable for persons with photosensitive conditions. Not suitable for light-sensitive epilepsy. Not suitable during pregnancy (unless approved by a doctor). If irritation, redness, or pain occurs, stop using the device immediately.

### Maintenance and Cleaning

Allow the device to cool down after use. Clean only with a soft dry or slightly damp cloth. Do not use alcohol, solvents, or aggressive cleaning agents. Do not immerse the device in water.

### Storage

Store in a dry place away from direct sunlight. Storage temperature: -10 °C to +45 °C. Protect from dust and humidity.

### Device Disposal

Do not dispose of with household waste. Take to a specialized collection point.

### Manufacturer

BlueTouch, s.r.o.  
Zelený pruh 1974/109,  
140 00 Prague, Czechia  
ID: 09437223, VAT: CZ09437223  
info@bluetouch.cz

## EU COMPLIANCE

This product has been designed and manufactured in accordance with the applicable legal regulations of the European Union relating to consumer electronics.

The product complies with the requirements of the following EU directives and regulations:

- Directive 2014/30/EU on electromagnetic compatibility (EMC),
- Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS),
- ČSN ISO 13485,
- ČSN EN IEC 55014-1.

Based on compliance with the above requirements, the product is marked with the CE marking.

The EU Declaration of Conformity is available from the importer and can be provided upon request to the relevant supervisory authorities or customers.

This product is intended exclusively for wellness and relaxation use.

It is not a medical device within the meaning of Regulation (EU) 2017/745 and is not intended for the diagnosis, treatment, or prevention of diseases.

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**MORE INFORMATION ABOUT LIGHT THERAPY  
CAN BE FOUND ON OUR WEBSITE.**

**BLUETOUCH.**